

PRODUCT BULLETIN


Product(s):	Fabius GS
Bulletin Number:	156
Release Date:	October 27, 2003
Pages:	1

Subject: Modification to Ventilator Hose-Moisture Reduction Kit

This product bulletin is a supplement to the previously released product bulletin #148 entitled "Ventilator Hose-Moisture Reduction Kit for Fabius GS."

General Overview

In order to provide additional flexibility in positioning the height of the watertrap, two different hose lengths are now included as part of the Moisture Reduction Kit. One hose, either 25cm or 35cm, can be used in position 2 shown in the diagram below. Selection of the hose length will depend upon the desired height at which the watertrap is to be located. To ensure optimal effectiveness, the watertrap should be located at the lowest point between the vent hoses. The addition of the second hose length should reduce the need for repositioning the Coby height to optimize the watertrap height.

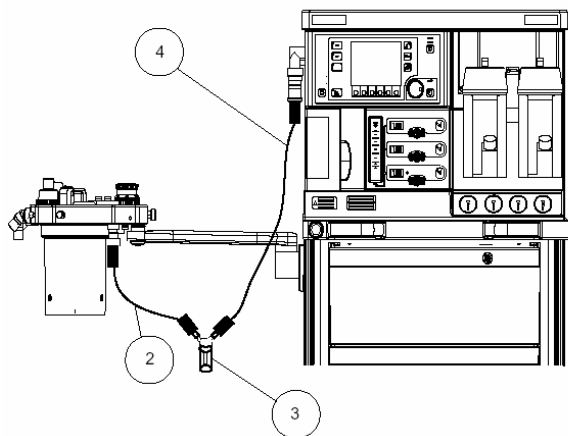

The part number and pricing of the moisture reduction kit remains the same.

<u>Description</u>	<u>Part Number</u>	<u>List Price</u>
Fabius GS Moisture Reduction Kit	4118299	\$150

Feel free to contact Mary Kyryk at ext. 2619 or kyrykm@draegermed.com if you have any additional questions.